


GLUTEAL FAT TRANSFER SET

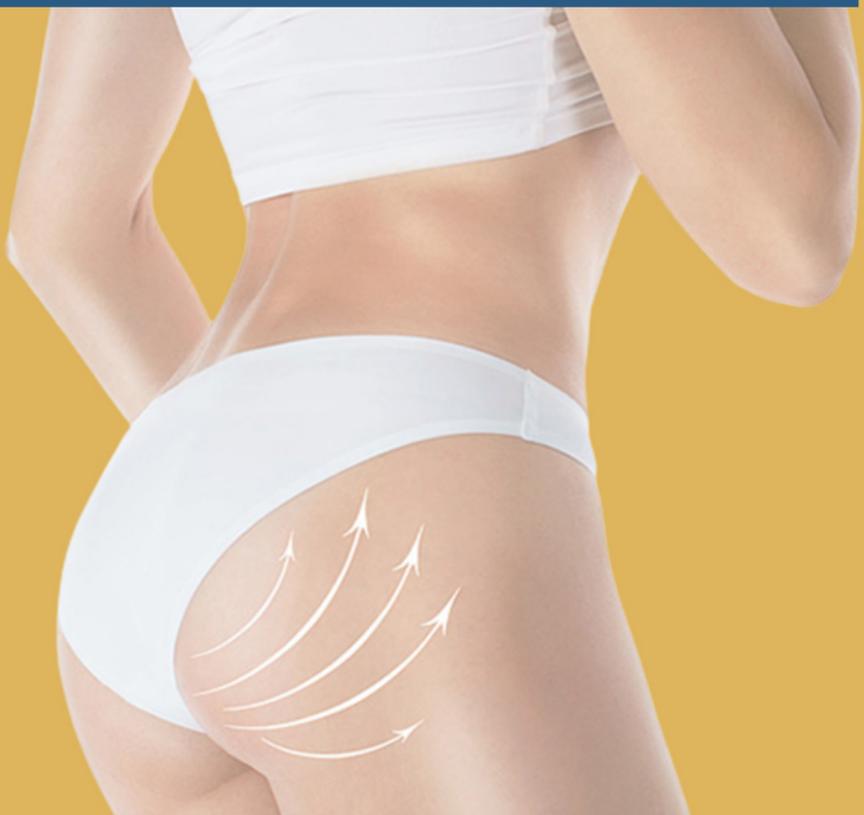

Redefining the Art of Fat Grafting

Gluteal Fat Transfer Set

Alira's Gluteal Fat Transfer Set containing reusable & cell-friendly cannulas designed by recommendation of experts perfect for fine gluteal sculpting and contouring. Alira's Cannulas are designed with your patient's safety in mind, so large volume fat grafting to the buttocks can be performed with confidence. This versatile set is perfect for Mommy Makeovers, gluteal augmentation, reshaping the mid-torso, and producing a unified, curvaceous shape of buttock.

1 Mercedes Liposuction Cannula - Three Lateral Ports

Alira's Mercedes Liposuction Cannula is a unique tool for liposuction procedures or fat harvesting. The three lateral ports aid to expedite the suction process in a variety of orientations. This circumferential three port configuration cell friendly liposuction cannula is effective for multi-plane soft tissue harvesting. In addition, this cannula is available with different sizes and working



Length: 20cm | **Diameter:** 3.5mm | **Qty:** 01

Length: 25cm | **Diameter:** 3.5mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub | Power Cannula | Fixed Handle | PAL Handle

2 Spoon Tip Injector-Blunt Tip

Alira's Spoon tip single port injector cannula with blunt tip design that passes through soft tissue effortlessly, while maintaining the integrity of the vascular structure. A spoon-shaped tip designed with a larger port for rapid injection of larger volumes, a groundbreaking breast and buttock workhorse. In addition, this cannula is available in different sizes and working lengths de-



Length: 15cm | **Diameter:** 3mm | **Qty:** 01

Length: 25cm | **Diameter:** 3mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub | Power Cannula | Fixed Handle | PAL Handle



3 Tumescent Infiltrator Straight

The safe and effective administration of tumescent solutions using tumescent infiltration cannulas ensures a painless liposuction procedure.



Length: 15cm | **Diameter:** 2.5mm (12G) | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub | Power Cannula | Fixed Handle | PAL Handle

4 Anaerobic Transfers

Alira's Anaerobic Transfers utilizing Micro Injector Cannulas are used to prepare and transport harvested tissue for fat grafting.



Type: Luer to Luer 2.4mm | **Qty:** 02

Type: Luer to Luer 1.4mm | **Qty:** 01

Type: Luer to Luer 1.2mm | **Qty:** 01

5 Decant Stands

To prepare and transfer harvested tissue for reinjection, anaerobic transfers are performed. Alira's Decant Stand's design also enables the decant of extracted tissue using natural gravity.

Type: 60cc to Luer | **Qty:** 01

Type: Luer Lock to Luer Lock | **Qty:** 01



6 Mesh For Nano Transfer

For separating big adipose tissue from small re-graftable fat cells through a Nano filter, Nano fat transfer mesh is excellent.



Size: 300 Micron | Qty: 01

Size: 500 Micron | Qty: 01

Size: 800 Micron | Qty: 01

7 Nano Transfer

Alira's Nano Transfer is a tool made to uniformly enlarge collected adipose tissue so that microneedles can readily inject it. The Nano Transfer has a cartridge. When doing liposuction and fat transfer operations, it is utilized to filter fat.

Qty: 01



WHY CHOOSE ALIRA?

Enhanced Precision: You can achieve precise fat graft placement, resulting in consistent and desirable cosmetic outcomes.

Minimized Patient Trauma: Alira's instruments minimize tissue trauma, leading to faster healing and reduced post-operative discomfort for patients.

Reduced Surgeon Fatigue: Ergonomic handles improve your comfort during procedures, reducing the risk of fatigue and enhancing surgical control.

Biocompatible Materials: All materials used are biocompatible, ensuring compatibility with the human body and reducing the risk of adverse reactions.



CORE AREAS OF EXPERTISE

Manufacturing
Extensive Range of
High Quality
Surgical Instruments

Instruments
Customization

Instruments
Repairing Services

Contract
Manufacturing

CONTACT US

ALIRA MEDICAL DEVICES LLC

+1 (818) 614 4592

info@alirameditives.com

www.alirameditives.com

14622 Ventura Blvd Ste 102, #712 Sherman Oaks, CA 91403

